



MEDITECH FILTERS:

MediStat™

The media used in the MediStat™ filter is electrostatically charged. This is a finely meshed membrane, which filters particles down to 2 microns in diameter.

The fibres in the filter are electrostatically charged with both positive and negative charge. This causes bacteria smaller than 2 microns to be attracted to the fibres.

All bacteria will be either positively or negatively charged - this results in high filtration without high resistance to flow, a major advantage over mechanical pleated filters.

All our filters reach 99.99% efficiency or greater

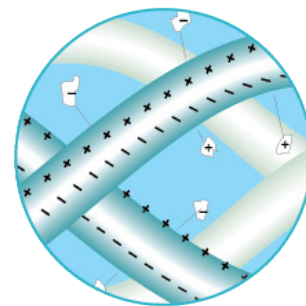


Illustration of the electrostatic process used in the MediStat™ filtration



Combination Filter Kits

- » 3 in 1 Combination Kits include a Filter, Catheter Mount and Elbow
- » 2 in 1 Kits include a Filter with a plain Elbow

Both kits assist clinical staff with set-ups for standard gas induction with all components provided in the one package

These components are all supplied with hard connections that make attachment and detachment simple and secure when in use.



222835XMMS



222836XMMS



222836PE

Product	Part No.	Box Quantity
MediStat™ Bacterial/Viral Filter, NO Gas Sample port, NO HME	222832	50
MediStat™ Bacterial/Viral Filter with Gas Sample port, NO HME	222832G	50
MediStat™ Bacterial/Viral Filter, with Gas Sample port, NO HME, Sterile	SU222832G	50
MediStat™ Bacterial/Viral Filter, with Gas Sample port, HME (Foam)	222835	50
MediStat™ Bacterial/Viral Filter, with Gas Sampling port, HME (Foam), Sterile	SU222835	50
MediStat™ Bacterial/Viral Filter, with Gas Sample port, HME (Paper)	222836	50
MediStat™ Bacterial/Viral Filter, with Gas Sample port, HME (Paper), Sterile	SU222836	50
MediStat™ Mini Bacterial/Viral Filter with Gas Sample port, HME (Paper)	222833	50
MediStat™ Mini Bacterial/Viral Filter with Gas Sample port, HME (Paper), Sterile	SU222833	50
MediStat™ Neo Bacterial/Viral Filter with Gas Sample port, with HME (Foam)	222832NEO	50
MediStat™ Neo Bacterial/Viral filter with Gas Sample port, with HME (Foam), Sterile	SU222832NEO	50
3 in 1 kit: Medistat™ 222835 Filter, Expandable Catheter Mount, Standard Elbow	222835XMMS	50
3 in 1 Kit: Medistat™ 222836 Filter, Expandable Catheter Mount, Standard Elbow	222836XMMS	50
2 in 1 Kit: Medistat™ 222836 Filter, Standard Elbow	222836PE	50

Contact our team for further information and support

JOHN PAXMAN
+61 (0) 434 361 558

MICHAEL SCHIRMER
+61 (0) 422 458 810

CUSTOMER SERVICE
+61 (0) 8 6162 9911