



IV INFUSION SET “ALL-IN-ONE”

Meditech's All-In-One IV Infusion Set is a high-quality, gravity fed IV line for everyday clinical use, that incorporates the functionality of a “chooks foot”.

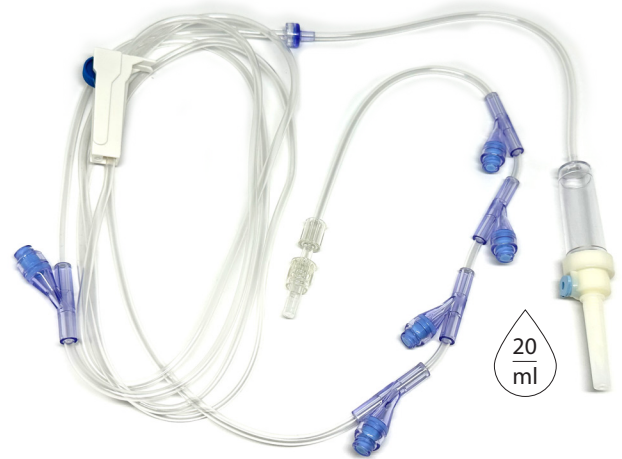

Features:

- » This particular IV set is 265cm in length. It incorporates a vented drip chamber with micron filter¹, a non-return valve², a roller clamp with IV tube holder³, 5 x needleless valve ports for injection⁴, and a rotating male luer at the distal tip⁵.
- » Incorporates 4 x lateral-laying, distal needless ports⁶ for multiple infusions, as required.
- » All-in-One: No leaks. No disconnection.
- » The IV set is STERILE, LATEX-free and DEHP-free. It provides 20 drops/ml.

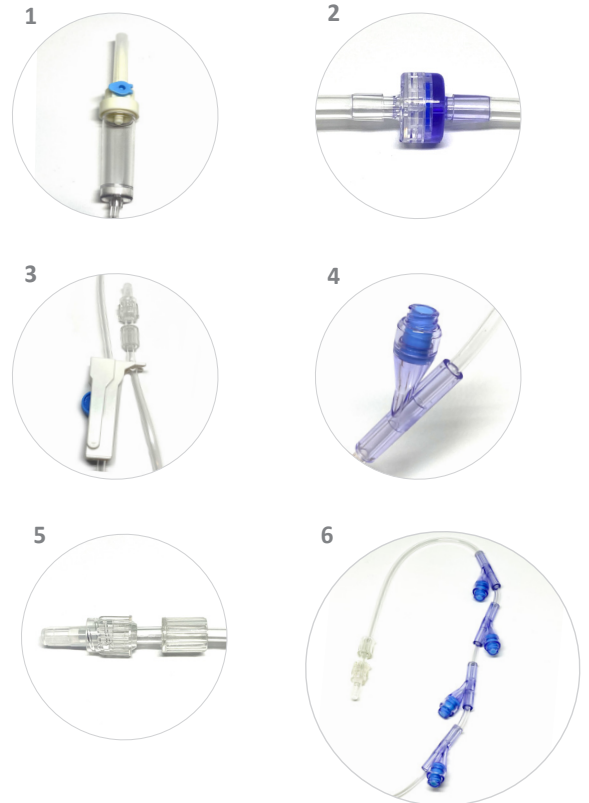

Product	Code	Box Quantity
IV Gravity All-In-One Infusion Set, 265cm, 5 x Needle-free Ports (4 x Distal Ports), Spiked Drip Chamber (20 drops/ml), Roller Clamp, Gravity-type backcheck valve.	88IV265B	300

Contact our team for further information and support

JOHN PAXMAN
+61 (0) 434 361 558

MICHAEL SCHIRMER
+61 (0) 422 458 810

CUSTOMER SERVICE
+61 (0) 8 6162 9911