

ANALGESIC BREATHING CIRCUITS

Meditech Systems Asia Pacific brings you a range of Entonox delivery circuits to meet the demands of pain management in Maternity and Emergency Departments, and other areas of the Hospital environment where Entonox gases are utilised.

Ergonomic design

Designed to allow the circuit to passively drape over the patients hand during administration, minimising the risk of the circuit being crushed or dropped.

One-Way Valves



Unique one-way valve to the inspiratory and expiratory ports of the circuits removes the need for an additional filter. This reduces resistance and minimises the required respiration effort of the patient.

Pink Scavenge with Connector

Pink tubing delineates the scavenge limb from the inspiratory limb. Includes a 19mm adaptor for easy connection to wall or mobile scavenge systems.

2.8m Expandable Tubing/Short Tubing

The long length of the inspiratory (and expiratory) limb permits patient movement, while other circuits allow the entonox system to be placed close to the administration site when required.

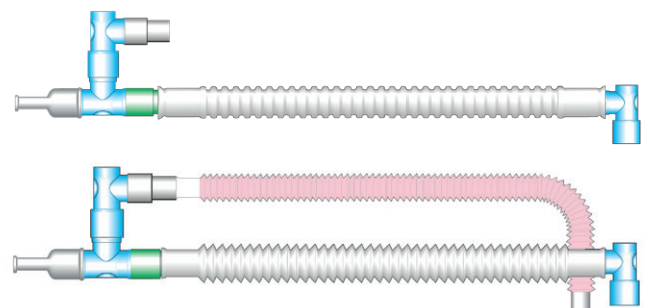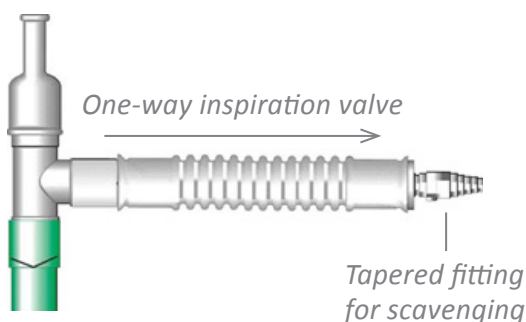

Product	Part No.	Box Quantity
Entonox Circuit, 2.8m, expandable, no Scavenge	22ENTX282200	Box of 20
Entonox Circuit, 2.8m, with mouthpiece + 2.8m pink Scavenge Limb tubing, 19mm Scavenge connector	22ENTX282228	Box of 20
Entonox unit with mouthpiece, one way inspiratory valve, and a 15cm tube	22ENV1508	Box of 20

Contact our team for further information and support

JOHN PAXMAN
+61 (0) 434 361 558

MICHAEL SCHIRMER
+61 (0) 422 458 810

CUSTOMER SERVICE
+61 (0) 8 6162 9911